Instructions for Use

TBE IgG ELISA

(Tick-borne Encephalitis Virus)

IVD

REF EIA-3860

Σ 96

DRG 

DRG Instruments GmbH, Germany
Frauenbergstraße 18, 35039 Marburg
Phone: +49 (0)6421-1700 0, Fax: +49 (0)6421-1700 50
Website: www.drg-diagnostics.de
E-mail: drg@drg-diagnostics.de

DRG 

DRG International, Inc., USA
841 Mountain Ave., Springfield, NJ 07081
Phone: (973) 564-7555, Fax: (973) 564-7556
Website: www.drg-international.com
E-mail: corp@drg-international.com

**Please use only the valid version of the Instructions for Use provided with the kit.
Verwenden Sie nur die jeweils gültige, im Testkit enthaltene, Gebrauchsanweisung.
Si prega di usare la versione valida delle istruzioni per l'uso a disposizione con il kit.
Por favor, se usa solo la version valida de la metodico técnico incluido aqui en el kit.
Utilisez seulement la version valide des Instructions d'utilisation fournies avec le kit.**

Table of Contents / Inhaltsverzeichnis / Tabella die Contenuti / Tabla de Contenidos / Sommaire

1	INTENDED USE.....	3
2	PRINCIPLE OF THE ASSAY	3
3	MATERIALS.....	3
4	STABILITY AND STORAGE	4
5	REAGENT PREPARATION	4
6	SAMPLE COLLECTION AND PREPARATION.....	4
7	ASSAY PROCEDURE	5
8	RESULTS.....	6
9	SPECIFIC PERFORMANCE CHARACTERISTICS	7
10	LIMITATIONS OF THE PROCEDURE.....	7
11	PRECAUTIONS AND WARNINGS	8
1	VERWENDUNGSZWECK.....	9
2	TESTPRINZIP	9
3	MATERIALIEN	9
4	STABILITÄT UND LAGERUNG	10
5	VORBEREITUNG DER REAGENZIEN.....	10
6	ENTNAHME UND VORBEREITUNG DER PROBEN	10
7	TESTDURCHFÜHRUNG	11
8	BERECHNUNG DER ERGEBNISSE	12
9	TESTMERKMALE	13
10	GRENZEN DES VERFAHRENS.....	13
11	SICHERHEITSMASSNAHMEN UND WARNHINWEISE	14
1	USO PREVISTO	16
2	PRINCIPIO DEL TEST.....	16
3	MATERIALI	16
4	MODALITÀ DI CONSERVAZIONE	17
5	PREPARAZIONE DEI REAGENTI.....	17
6	PRELIEVO E PREPARAZIONE DEI CAMPIONI	17
7	PROCEDIMENTO	18
8	RISULTATI.....	19
9	CARATTERISTICHE DEL TEST.....	20
10	LIMITAZIONI	20
11	PRECAUZIONI E AVVERTENZE.....	21

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1 INTENDED USE

The TBE IgG ELISA is intended for the quantitative determination of IgG class antibodies against TBE in human serum or plasma (citrate, heparin).

2 PRINCIPLE OF THE ASSAY

The quantitative immunoenzymatic determination of specific antibodies is based on the ELISA (Enzyme-linked Immunosorbent Assay) technique.

Microtiterplates are coated with specific antigens to bind corresponding antibodies of the sample. After washing the wells to remove all unbound sample material a horseradish peroxidase (HRP) labelled conjugate is added. This conjugate binds to the captured antibodies. In a second washing step unbound conjugate is removed. The immune complex formed by the bound conjugate is visualized by adding Tetramethylbenzidine (TMB) substrate which gives a blue reaction product.

The intensity of this product is proportional to the amount of specific antibodies in the sample. Sulphuric acid is added to stop the reaction. This produces a yellow endpoint colour. Absorbance at 450/620 nm is read using an ELISA Microtiterplate reader.

3 MATERIALS

3.1 Reagents supplied

- **Microtiterplate (IgG):**
12 break apart 8-well snap-off strips coated with TBE/FSME Virus antigens; in resealable aluminium foil.
- **DIL:**
1 bottle containing 100 mL of phosphate buffer (10 mM) for sample dilution; pH 7.2 ± 0.2; coloured yellow; ready to use; white cap; ≤ 0.0015% (v/v) CMIT/ MIT (3:1).
- **SOLN STOP:**
1 bottle containing 15 mL sulphuric acid, 0.2 mol/L; ready to use; red cap.
- **WASH BUF 20x:**
1 bottle containing 50 mL of a 20-fold concentrated phosphate buffer (0.2 M), pH 7.2 ± 0.2, for washing the wells; white cap; 0.2% (w/v) 5-Bromo-5-nitro-1,3-dioxane.
- **Conjugate:**
1 bottle containing 20 mL of peroxidase labelled antibody to human IgG in phosphate buffer (10 mM); coloured blue; ready to use; black cap.
- **SUB TMB:**
1 bottle containing 15 mL 3,3',5,5'-tetramethylbenzidine (TMB), < 0.1%; ready to use; yellow cap.
- **Control:**
1 bottle containing 2 mL; coloured yellow; ready to use; orange cap; ≤ 0.02% (v/v) MIT.
- **Standards:**
5 vials, each containing 2 mL standard (human serum or plasma); coloured yellow; ready to use. ; ≤ 0.02% (v/v) MIT
Standard A: 0 DU/mL; blue cap
Standard B: 50 DU/mL; green cap
Standard C: 130 DU/mL; yellow cap [DU = DRG Units]
Standard D: 200 DU/mL; red cap
Standard E: 300 DU/mL; white cap

For hazard and precautionary statements see 11.1

3.2 Materials supplied

- 1 Cover foil
- 1 Instructions for use (IFU)

3.3 Materials and Equipment needed

- ELISA Microtiterplate reader, equipped for the measurement of absorbance at 450/620 nm
- Incubator 37 °C
- Manual or automatic equipment for rinsing Microtiterplates
- Pipettes to deliver volumes between 10 and 1000 µL
- Vortex tube mixer
- Distilled water
- Disposable tubes

4 STABILITY AND STORAGE

Store the kit at 2 °C - 8 °C.

The opened reagents are stable up to the expiry date stated on the label when stored at 2 °C - 8 °C.

5 REAGENT PREPARATION

It is very important to bring all reagents and samples to room temperature (20 °C - 25 °C) and mix them before starting the test run!

5.1 Microtiterplate

The break-apart snap-off strips are coated with TBE/FSME Virus antigens.

Immediately after removal of the strips, the remaining strips should be resealed in the aluminium foil along with the desiccant supplied and stored at 2 °C - 8 °C.

5.2 WASH BUF 20x

Dilute WASH BUF 20x 1 + 19; e. g. 10 mL WASH BUF 20x + 190 mL distilled water.

The diluted buffer (WASH BUF 1x) is stable for 5 days at room temperature (20 °C - 25 °C). In case crystals appear in the concentrate, warm up the solution to 37 °C e.g. in a water bath. Mix well before dilution.

5.3 SUB TMB

The reagent is ready to use and has to be stored at 2 °C - 8 °C, away from the light.

SUB TMB should be colourless or could have a slight blue tinge. If SUB TMB turns into blue, it may have become contaminated and should be thrown away.

6 SAMPLE COLLECTION AND PREPARATION

Use human serum or plasma (citrate, heparin) samples with this assay.

If the assay is performed within 5 days after sample collection, the samples should be kept at 2 °C - 8 °C; otherwise they should be aliquoted and stored deep-frozen (-70 °C to -20 °C). If samples are stored frozen, mix thawed samples well before testing. Avoid repeated freezing and thawing.

Heat inactivation of samples is not recommended.

6.1 Sample Dilution

Before assaying, all samples should be diluted 1 + 100 with DIL.

Dispense 10 µL sample and 1 mL DIL into tubes to obtain a 1+100 dilution and thoroughly mix with a Vortex.

7 ASSAY PROCEDURE

Please read the instructions for use carefully **before** performing the assay. Result reliability depends on strict adherence to the instructions for use as described. The following test procedure is only validated for manual procedure. If performing the test on ELISA automatic systems we recommend increasing the washing steps from three to five and the volume of **WASH|BUF|1x** from 300 µL to 350 µL to avoid washing effects. Pay attention to chapter 11. Prior to commencing the assay, the distribution and identification plan for all samples and standards/controls (duplicates recommended) should be carefully established. Select the required number of microtiter strips or wells and insert them into the holder.

Perform all assay steps in the order given and without any delays.

A clean, disposable tip should be used for dispensing each standard/control and sample.

Adjust the incubator to 37 °C ± 1 °C.

1. Dispense 100 µL standards/controls and diluted samples into their respective wells. Leave well A1 for the Substrate Blank.
2. Cover wells with the foil supplied in the kit.
3. **Incubate for 1 hour ± 5 min at 37 °C ± 1 °C.**
4. When incubation has been completed, remove the foil, aspirate the content of the wells and wash each well three times with 300 µL of **WASH|BUF|1x**. Avoid overflows from the reaction wells. The interval between washing and aspiration should be > 5 sec. At the end carefully remove remaining fluid by tapping strips on tissue paper prior to the next step!
Note: Washing is important! Insufficient washing results in poor precision and false results.
5. Dispense 100 µL Conjugate into all wells except for the Substrate Blank well A1.
6. **Incubate for 30 min at room temperature (20 °C - 25 °C).** Do not expose to direct sunlight.
7. Repeat step 4.
8. Dispense 100 µL **SUB|TMB** into all wells.
9. **Incubate for exactly 15 min at room temperature (20 °C - 25 °C) in the dark.** A blue colour occurs due to an enzymatic reaction.
10. Dispense 100 µL **SOLN|STOP** into all wells in the same order and at the same rate as for **SUB|TMB**, thereby a colour change from blue to yellow occurs.
11. Measure the absorbance at 450/620 nm within 30 min after addition of the **SOLN|STOP**.

7.1 Measurement

Adjust the ELISA Microtiterplate reader **to zero** using the **Substrate Blank**.

If - due to technical reasons - the ELISA Microtiterplate reader cannot be adjusted to zero using the Substrate Blank, subtract its absorbance value from all other absorbance values measured in order to obtain reliable results!

Measure the absorbance of all wells at **450 nm** and record the absorbance values for each standard/control and sample.

Bichromatic measurement using a reference wavelength of 620 nm is recommended.

Where applicable calculate the **mean absorbance values** of all duplicates.

8 RESULTS

8.1 Run Validation Criteria

In order for an assay run to be considered valid, these instructions for use have to be strictly followed and the following criteria must be met:

- **Substrate blank:** Absorbance value < 0.100
- **Standard A:** Absorbance value < 0.200
- **Standard B:** Absorbance value > 0.050
- **Standard C:** Absorbance value > Standard B
- **Standard D:** Absorbance value > Standard C
- **Standard E:** Absorbance > 1.000
- **Control:** Result in DU/mL within range indicated on the label

Standard A < Standard B < Standard C < Standard D < Standard E

If these criteria are not met, the test is not valid and must be repeated.

8.2 Calculation of Results

In order to obtain quantitative results in DU/mL plot the (mean) absorbance values of the 5 Standards A, B, C, D and E on (linear/Linear) graph paper in a system of coordinates against their corresponding concentrations (0, 50, 130, 200 and 300 DU/mL) and draw a standard calibration curve (absorbance values on the vertical y-axis, concentrations on the horizontal x-axis).

Read results from this standard curve employing the (mean) absorbance values of each patient sample and control.

For the calculation of the standard-curve mathematical Point to Point function should be used.

8.3 Interpretation of Results

Normal value ranges for this ELISA should be established by each laboratory based on its own patient populations in the geographical areas serviced.

The following values should be considered as a guideline:

Positive	> 110 DU/mL	Antibodies against the pathogen are present. There has been a contact with the antigen (pathogen resp. vaccine). <u>After vaccination:</u> This is a case of seroconversion. Check anamnestic data and if necessary complete basic immunization or give a booster
Equivocal	55 –110 DU/mL	Antibodies against the pathogen could not be detected clearly. It is recommended to repeat the test with a fresh sample in 2 to 4 weeks. <u>After vaccination:</u> This may be a case of seroconversion. Continue with basic immunisation or booster. A non-specific reaction is not excluded.
Negative	< 55 DU/mL	The sample contains no antibodies against the pathogen. A previous contact with the antigen (pathogen resp. vaccine) is unlikely. <u>After vaccination:</u> No seroconversion after vaccination. This may be the case after the first vaccination during basic immunisation. Patients which show low or no response.
Diagnosis of an infectious disease should not be established on the basis of a single test result. A precise diagnosis should take into consideration clinical history, symptomatology as well as serological data. In immunocompromised patients and newborns serological data only have restricted value.		

8.3.1 Antibody Isotypes and State of Infection

Serology	Significance
IgM	Characteristic of the primary antibody response High IgM titer with low IgG titer: → suggests a current or very recent infection Rare: → persisting IgM
IgG	Characteristic of the secondary antibody response May persist for several years High IgG titer with low IgM titer: → may indicate a past infection

9 SPECIFIC PERFORMANCE CHARACTERISTICS

The results refer to the groups of samples investigated; these are not guaranteed specifications.

9.1 Precision

Intra assay	n	Mean (E)	CV (%)
#1	24	1.411	3.51
#2	24	1.991	2.23
#3	24	1.913	8.93
Inter assay	n	Mean (DU/mL)	CV (%)
#1	12	26.55	9.18
#2	12	15.38	12.99
#3	12	253.81	10.48

9.2 Diagnostic Specificity

The diagnostic specificity is defined as the probability of the assay of scoring negative in the absence of the specific analyte.

It is 100% (95% confidence interval: 93.84% - 100%).

9.3 Diagnostic Sensitivity

The diagnostic sensitivity is defined as the probability of the assay of scoring positive in the presence of the specific analyte.

It is 100% (95% confidence interval: 96.31% - 100%).

9.4 Analytical Sensitivity

The analytical sensitivity (according to CLSI EP17-A) is defined as the apparent concentration of the analyte that can be distinguished from the zero calibrator. It is 1.78 DU/mL.

9.5 Interferences

Interferences with hemolytic, lipemic or icteric samples are not observed up to a concentration of 10 mg/mL hemoglobin, 5 mg/mL triglycerides and 0.5 mg/mL bilirubin.

9.6 Cross Reactivity

Cross reactivity with other flaviviruses cannot be excluded and should be taken into account for result interpretation.

9.7 Measurement range

The measurement range is 1.78 DU/mL – 300 DU/mL.

10 LIMITATIONS OF THE PROCEDURE

Bacterial contamination or repeated freeze-thaw cycles of the sample may affect the absorbance values.

11 PRECAUTIONS AND WARNINGS

- The test procedure, the information, the precautions and warnings in the instructions for use have to be strictly followed. The use of the test kits with analyzers and similar equipment has to be validated. Any change in design, composition and test procedure as well as for any use in combination with other products not approved by the manufacturer is not authorized; the user himself is responsible for such changes. The manufacturer is not liable for false results and incidents for these reasons. The manufacturer is not liable for any results by visual analysis of the patient samples.
- Only for in-vitro diagnostic use.
- All materials of human or animal origin should be regarded and handled as potentially infectious.
- All components of human origin used for the production of these reagents have been tested for anti-HIV antibodies, anti-HCV antibodies and HBsAg and have been found to be non-reactive.
- Do not interchange reagents or Microtiterplates of different production lots.
- No reagents of other manufacturers should be used along with reagents of this test kit.
- Do not use reagents after expiry date stated on the label.
- Use only clean pipette tips, dispensers, and lab ware.
- Do not interchange screw caps of reagent vials to avoid cross-contamination.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- After first opening and subsequent storage check conjugate and standard/control vials for microbial contamination prior to further use.
- To avoid cross-contamination and falsely elevated results pipette patient samples and dispense reagents without splashing accurately into the wells.
- The ELISA is only designed for qualified personnel following the standards of good laboratory practice (GLP).
- For further internal quality control each laboratory should additionally use known samples.

11.1 Safety note for reagents containing hazardous substances

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Therefore, the following hazard and precautionary statements apply.

Warning

H317	May cause an allergic skin reaction.
P261	Avoid breathing spray
P280	Wear protective gloves/ protective clothing.
P302+P352	IF ON SKIN: Wash with plenty of soap and water.
P333+P313	If skin irritation or rash occurs: Get medical advice/ attention.
P362+P364	Take off contaminated cloths and wash it before reuse.

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Therefore, the following hazard and precautionary statements apply.

Warning

H315	Causes skin irritation.
H319	Causes serious eye irritation
P280	Wear protective gloves/ protective clothing.
P302+P352	IF ON SKIN: Wash with plenty of soap and water.
P305+P351+P338	IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.
P337+P313	If eye irritation persists: Get medical advice/attention.

Further information can be found in the safety data sheet.

11.2 Disposal Considerations

Residues of chemicals and preparations are generally considered as hazardous waste. The disposal of this kind of waste is regulated through national and regional laws and regulations. Contact your local authorities or waste management companies which will give advice on how to dispose hazardous waste.

For information about the packaging materials refer to PACKAGING MATERIALS.

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13 SCHEME OF THE ASSAY

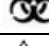
Test Preparation

<p>Prepare reagents and samples as described.</p> <p>Establish the distribution and identification plan for all samples and standards/controls.</p> <p>Select the required number of microtiter strips or wells and insert them into the holder.</p>
--

Assay Procedure

	Substrate Blank (A1)	Std A	Std B	Std C	Std D	Std E	Control	Sample (diluted 1+100)
Standard A	-	100 µL	-	-	-	-	-	-
Standard B	-	-	100 µL	-	-	-	-	-
Standard C	-	-	-	100 µL	-	-	-	-
Standard D	-	-	-	-	100 µL	-	-	-
Standard E	-	-	-	-	-	100 µL	-	-
Control	-	-	-	-	-	-	100 µL	-
Sample (diluted 1+100)	-	-	-	-	-	-	-	100 µL
Cover wells with foil supplied in the kit Incubate for 1 h at 37 °C ± 1 °C Wash each well three times with 300µl of WASH BUF 1x								
Conjugate	-	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
Incubate for 30 min at room temperature (20 °C - 25 °C) Do not expose to direct sunlight Wash each well three times with 300 µL of WASH BUF 1x								
SUB TMB	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
Incubate for exactly 15 min at room temperature (20 °C - 25 °C) in the dark								
SOLN STOP	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL	100 µL
Photometric measurement at 450 nm (reference wavelength: 620 nm)								

SYMBOLS USED

Symbol	English	Deutsch	Italiano	Español	Français
	European Conformity	CE-Konformitätskennzeichnung	Conformità europea	Conformidad europea	Conformité normes européennes
	Consult instructions for use *	Gebrauchsanweisung beachten *	Consultare le istruzioni per l'uso	Consulte las instrucciones de uso	Consulter les instructions d'utilisation
IVD	<i>In vitro</i> diagnostic medical device *	<i>In-vitro</i> -Diagnostikum *	Dispositivo medico-diagnostico in vitro	Producto sanitario para diagnóstico In vitro	Dispositif médical de diagnostic in vitro
REF	Catalogue number *	Artikelnummer *	Numero di Catalogo	Número de catálogo	Référence de catalogue
LOT	Batch code *	Fertigungslosnummer, Charge *	Codice del lotto	Código de lote	Numéro de lot
	Contains sufficient for <n> tests *	Ausreichend für <n> Prüfungen *	Contenuto sufficiente per "n" saggi	Contenido suficiente para <n> ensayos	Contenu suffisant pour "n" tests
	Temperature limit *	Temperaturbegrenzung *	Temperatura di conservazione	Temperatura de conservación	Température de conservation
	Use-by date *	Verwendbar bis *	Utilizzare prima del	Establa hasta	Utiliser jusque
	Manufacturer *	Hersteller *	Fabbricante	Fabricante	Fabricant
	Biological risks*	Biologische Risiken*	Rischi biologici	Riesgos biológicos	Risques biologiques
	Caution *	Achtung *	Attenzione	Precaución	Attention
UDI	Unique Device Identifier	Eindeutige Produktidentifizierung	identificazione unica del dispositivo	identificación única del producto	identification unique des dispositifs
RUO	For research use only	Nur für Forschungszwecke	Solo a scopo di ricerca	Sólo para uso en investigación	Seulement dans le cadre de recherches
	Distributed by	Vertreiber	Distributore	Distribuidor	Distributeur
<i>Content</i>	Content	Inhalt	Contenuto	Contenido	Contenu
<i>Volume/No.</i>	Volume / No.	Volumen/Anzahl	Volume/Quantità	Volumen/Número	Volume/Quantité
MTP	Microtiterplate	Mikrotiterplatte	Piastre di Microtitolazione	Placa de Microtitulación	Plaque de Microtitrage
CONJ	Conjugate	Konjugat	Coniugato	Conjugado	Conjugué
CONTROL -	Negative Control	Negativkontrolle	Controllo negativo	Control positivo	Contrôle négatif
CONTROL +	Positive Control	Positivkontrolle	Controllo positivo	Control positivo	Contrôle positif
CUT OFF	Cut off control	Cut off-Kontrolle	Controllo Cut-off	Control Cut-off	Contrôle Cut-off
CAL	Standard or Calibrator	Standard oder Kalibrator	Standard o Calibratore	Estándar o Calibrador	Standard o Etalon
DIL	Sample Dilution Buffer	Probenverdünnungspuffer	Tampone di Diluizione del Campione	Tampón de Dilución de Muestra	Tampon de Dilution d'Échantillon
DIL A	IgA Sample Dilution Buffer	IgA-Probenverdünnungspuffer	Tampone di Diluizione del Campione IgA	Tampón de Dilución de Muestra IgA	Tampon de Dilution d'Échantillon IgA
DIL M	IgM Sample Dilution Buffer	IgM-Probenverdünnungspuffer	Tampone di Diluizione del Campione IgM	Tampón de Dilución de Muestra IgM	Tampon de Dilution d'Échantillon IgM
SOLN STOP	Stop solution	StoppLösung	Soluzione bloccante	Solución de parada	Solution d'arrêt
SUB TMB	TMB Substrate solution	TMB-Substratlösung	Soluzione substrato TMB	Solución de sustrato de TMB	Solution de substrat TMB
WASH BUF 20x	Washing Buffer 20x concentrated	Waschpuffer 20x konzentriert	Tampone di lavaggio concentrazione x20	Tampón de lavado concentrado x20	Tampon de lavage concentré 20 x
WASH BUF 1x	20-fold dilution of WASH BUF 20x	20-fach Verdünnung von WASH BUF 20x	Diluizione 20 volte del WASH BUF 20x	Dilución de 20 veces del WASH BUF 20x	Dilution 20 fois du WASH BUF 20x